#407

PAVIS S.p.A.

Sede / Headquarter:
Via G. Rossini, 11/B - 21020 Buguggiate (VA) - Italy
Tel. (+39) 0332 45 54 63
Fax (+39) 0332 45 97 48

www.pavis.it - info@pavis.it
pavispec@certimprese.it
R.I. / C.F. / P.IVA 01286640121
R.N. / F.C. / VAT IT01286640121

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
(UK)

1 - INDICATIONS
- Traumatic pathologies: recurrent shoulder dislocation and subluxation
- Neuro-degenerative diseases: osteoarthritis, degenerative processes of the shoulder, neurological shoulder
- Post surgery: Outcomes of surgeries for glenohumeral joint instability, shoulder dislocation and shoulder arthroscopy
- As prevention of shoulder dislocations/subluxations or instability, both in sports and in daily practice, conservative treatment of the neurological shoulder

2 - CONTRA-INDICATIONS AND SIDE-EFFECTS
The choice between different models must be done according to the pathology and the necessities stated by the doctor.
Use should be limited to the pathology period. Do not use the device if it create pressure on wounds and swellings. Do not wear it close to open flames. To avoid skin redness (caused by movements mechanical frictions) it is preferable to wear it over a cotton t-shirt.
Since this device is NON-INVASIVE and made with non-allergenic fabrics, no negative or side effects have ever been noted.
Any serious incident occurred in relation to the device use, should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

3 - MAINTENANCE AND CARE
• Before washing, close hook and loop parts

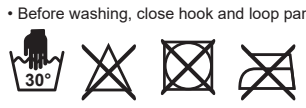This device is non-biodegradable. Even though it is made of natural fabrics, after use, avoid release the product to the environment. AT THE END OF ITS USABLE LIFE HAS TO BE DISCHARGED ALONG WITH SOLID ORDINARY WASTE.

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