#574-584

PAVIS S.p.A.

Sede / Headquarter:
Via G. Rossini, 11/B - 21020 Buguggiate (VA) - Italy
Tel. (+39) 0332 45 54 63
Fax (+39) 0332 45 97 48

www.pavis.it - info@pavis.it
pavispec@certimprese.it
R.I. / C.F. / P.IVA 01286640121
R.N. / F.C. / VAT IT01286640121

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
(UK)

1 - INDICATIONS
- Lumbago, ischialgia, lumbar protrusions, lumbar arthrosis, lumbar muscle contractures, discopathy, lumbar hyper/hypo-lordosis, osteoporosis, crural pains (with lumbar origin), vertebral collapse, spondyloarthrosis

2 - CONTRA-INDICATIONS AND SIDE-EFFECTS
The choice between different models must be done according to the pathology and the necessities stated by the doctor. To avoid weaken back muscles, it is recommended not to use the corset more than time needed to solve the pathology. Do not use the device if it create pressure on wounds and swellings. Do not wear it close to open flames.
Since this device is NON-INVASIVE and made with non-allergenic fabrics, no negative or side effects have ever been noted.
Any serious incident occurred in relation to the device use, should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

The first custom application of this device must be performed by a professional.
Failure to do so will void the warranty and reduce the benefical effects.
Remove the panel from its cover and custom mold it using either a hairdrier, an oven or by dipping it in hot water.
CAUTION: TO AVOID SERIOUS BURNS, MAKE SURE YOU DO NOT OVERHEAT THE PANEL.

3 - MAINTENANCE AND CARE
- Before washing, remove the panel - fold the Lumbosacral Brace and fasten the hook and loop parts

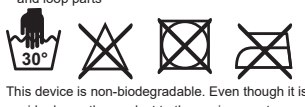This device is non-biodegradable. Even though it is made of natural fabrics, after use, avoid release the product to the environment.
AT THE END OF ITS USABLE LIFE HAS TO BE DISCHARGED ALONG WITH SOLID ORDINARY WASTE.

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